



HEMO-CATH® LT

Indications for Use: The Medcomp® Hemo-Cath® LT Silicone Double Lumen Catheter can be utilized for long term implantation as well as temporary access for hemodialysis, hemoperfusion, or apheresis therapy. The cannula may be inserted via the Seldinger Technique due to the inner Teflon stylet, increasing linear strength. The stylet is removed after insertion, leaving the soft silicone cannula in the body. The flexible silicone make-up conforms well to the vessel anatomy, resulting in higher patient tolerance during extended use.


Contraindications: The Subclavian Approach is NOT recommended for use with the Medcomp® Double Lumen Subclavian-Femoral Catheter in Hemodialysis or Hemoperfusion Procedures used for the management of acute poisoning or other situations in which a ventilator might be used due to risk of traumatic pneumothorax posing a dangerous complication for the patient.

Refer to Instructions for Use provided with the product for complete instructions, warnings, precautions, and contraindications. Observe all instructions for use prior to using products. Failure to do so may result in patient complications.



CIRCLE "C" INTERNAL LUMEN DESIGN

Silicone Material	Step Tip Tip Design	8F, 12.5F French Sizes	Straight, Pre-Curved Configurations	Antegrade Insertion Type
-----------------------------	-------------------------------	----------------------------------	---	------------------------------------

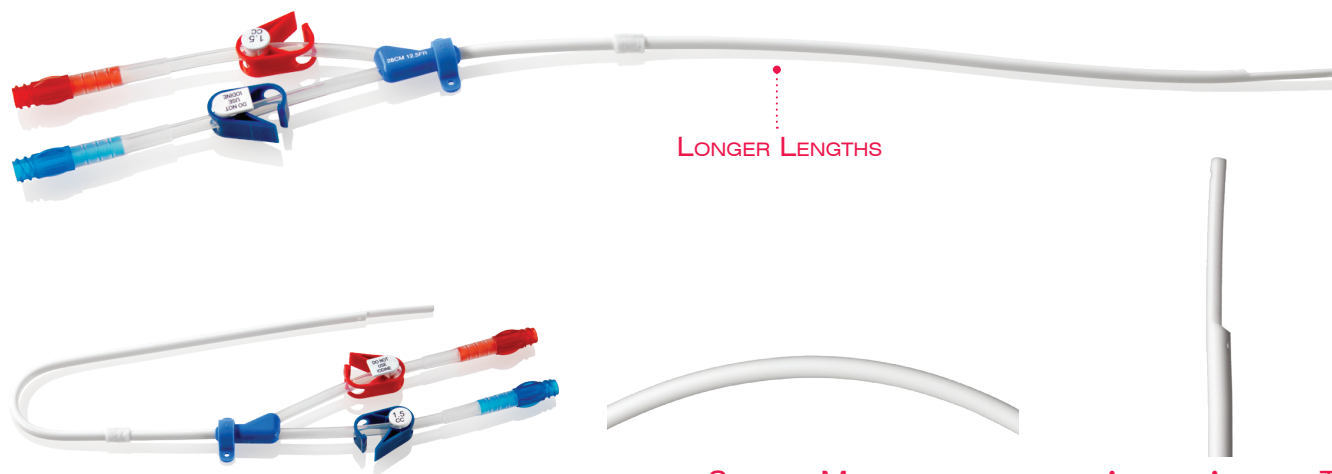

LONGER LENGTHS

SILICONE MATERIAL

ANGLED ARTERIAL TIP

Refer to the Table of Contents for Important Risk Information regarding this device.

PN2695 Medcomp Catalog DOM Rev. B 1/21 Excerpt

8F BASIC SET		
DSL18P	8F X 18CM HEMO-CATH® LT CATHETER, STRAIGHT	5/BOX
DSL24P	8F X 24CM HEMO-CATH® LT CATHETER, STRAIGHT	5/BOX

12.5F BASIC SET		
DMC101241	12.5F X 15CM HEMO-CATH® LT CATHETER, STRAIGHT	5/BOX
DMC101242	12.5F X 18CM HEMO-CATH® LT CATHETER, STRAIGHT	5/BOX
DMC101243	12.5F X 24CM HEMO-CATH® LT CATHETER, STRAIGHT	5/BOX
DSL28	12.5F X 28CM HEMO-CATH® LT CATHETER, STRAIGHT	5/BOX
DSL32	12.5F X 32CM HEMO-CATH® LT CATHETER, STRAIGHT	5/BOX
DSL28PC	12.5F X 28CM HEMO-CATH® LT CATHETER, PRE-CURVED	5/BOX
DSL32PC	12.5F X 32CM HEMO-CATH® LT CATHETER, PRE-CURVED	5/BOX

12.5F FULL TRAY		
MCDLTSL28	12.5F X 28CM HEMO-CATH® LT CATHETER, STRAIGHT	5/BOX
MCDLTSL32	12.5F X 32CM HEMO-CATH® LT CATHETER, STRAIGHT	5/BOX
MCDLTSL28PC	12.5F X 28CM HEMO-CATH® LT CATHETER, PRE-CURVED	5/BOX
MCDLTSL32PC	12.5F X 32CM HEMO-CATH® LT CATHETER, PRE-CURVED	5/BOX

8F BASIC SET CONTENTS:

- (1) Catheter
- (1) Introducer Needle
- (1) Guidewire
- (1) Tunneler
- (1) Peelable Introducer
- (1) Syringe
- (1) Hemo-Cath® Hub Clip
- (1) Adhesive Wound Dressing
- (2) End Caps
- (1) Patient ID Card

12.5F BASIC SET CONTENTS:

- (1) Catheter
- (1) Introducer Needle
- (1) Guidewire
- (1) Scalpel
- (1) Peelable Introducer
- (1) Tunneler
- (1) Hemo-Cath® Hub Clip
- (1) Adhesive Wound Dressing
- (2) End Caps
- (1) Patient ID Card

12.5F FULL TRAY CONTENTS:

- (1) Catheter
- (1) Introducer Needle
- (1) Guidewire
- (1) Tunneler
- (1) Peelable Introducer
- (2) 5cc Syringe
- (2) 10cc Syringe
- (1) 5cc Lidocaine
- (1) Scalpel
- (3) 4" x 4" Gauze
- (1) Surgical Gloves
- (1) Fenestrated Drape
- (1) Hemostat
- (3) Hypodermic Needles
- (1) Chloraprep
- (1) 2-0 Silk Suture
- (1) Hemo-Cath® Hub Clip
- (1) Adhesive Wound Dressing
- (2) End Caps

Refer to the Table of Contents for Important Risk Information regarding this device.